



**BANBROS ENGINEERING PVT. LTD.**  
The Precision Measurement People

# PRESSURE COLD MOUNTING MACHINE

- **VALUE**
- **VERSATILITY**
- **PERFORMANCE**



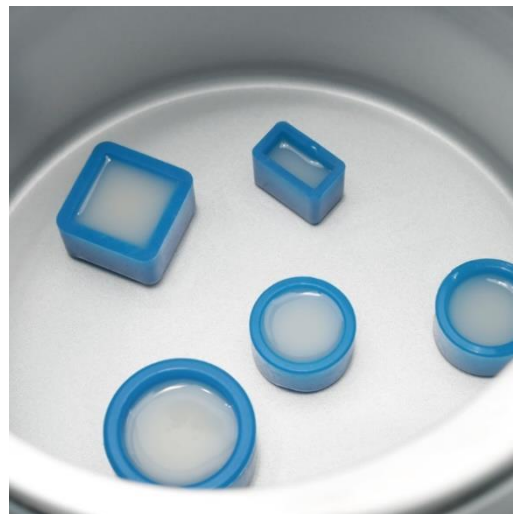
**BTM-P**

**DESCRIPTION:**

Our transparent moulding machine is suitable for the sample preparation of metallographic cold mounting, and better for epoxy or acrylic cold mounting.

**FEATURES:**

- No electricity, pressurized by a pneumatic system
- With water-separator and fast-insert connectors
- Quick lock and release cover design.
- High performance in sealing, and overflow valve in safety.
- More spacious mounting chamber, and stainless-steel container.
- Equipped with a timer to remind curing time

**TECHNICAL SPECIFICATIONS:**

Model No.	BTM-P
Operating Pressure	2.0 bar
Safety Valve Overflow Pressure	2.8 ~ 3.3. bar
Air Pressure Inlet Pressure	3~10 bar
Chamber Diameter	Φ235×170mm
Size (W x D x H) in mm	340 x 400 x 270 (Closed)
	340 x 400 x 370 (Open)
Weight	Approx. 8.5 Kg